



Elite Pro M810A

Q Switch ND YAG Tattoo & Pigmentation Removal Machine

Application





1064nm Tip (Standard)

Professional for dark colors tattoo, pigment removal, nail fungus.





532nm Tip (Standard)

Specialized in colorized tattoo treatment, such as red, brown and coffee etc.





1320nm Tip (Standard)

Mainly applied to deep cleaning, skin whitening, pore shrinking.






Fractional Tip (Optional)

Skin reconstruction, skin rejuvenation



The specific laser wavelength effect on the tattoo pigments and melanin of the epidermis and dermis. The melanin and pigments begin to swell and explode to tiny particles after absorbing this high energy laser. One part of them will be shot out of the body, the other parts will be gradually removed by metabolism.

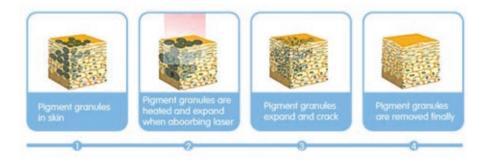

Proven Clinical Results



Advantages



1. TUV Medical CE approved

Strict quality testing process, safe treatment and stable working.



2. Portable but powerful performance

Double laser rods with super energy 1600mj, spot size 2 times larger, high efficiency and 50% treatment time saved.



3. High-efficiency Dermal Pigment & Tattoo Removal

Double energy handpiece especially focus on dermal pigment and colorful tattoo removal, permanent results and less sessions is needed.



4. Super long lifetime

Lifetime of xenon lamp more than 1,000,000 shots.



5. Intelligent interface

16 kinds of languages available for more convenient operation.

■ Technique Specification

Model NO.	Med-810A
wavelength	1064nm, 532nm, 1320nm standard;
	Fractional tip optional
Laser Power	500W
Energy	300-1600mj
Yag bar size	Φ7+φ8
Display	8.4" True Color LCD Touch Screen
Cooling	Air+Water+Twin Radiators Cooling System
Net Weight	25kg
Dimensions (WxLxH)	45*40*57cm